FCC Caution.

§ 15.19 Labelling requirements.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

§ 15.21 Information to user.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

§ 15.105 Information to the user.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

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RSS-GEN Issue 4 November 2014 "&" Rss-gen délivré le 4 novembre 2014

- English:

This device complies with Industry Canada licence-exempt RSS standard(s).

Operation is subject to the following two conditions: (1) This device may not cause interference, and (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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
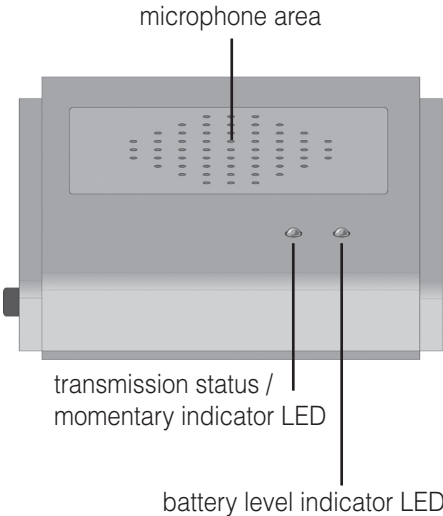
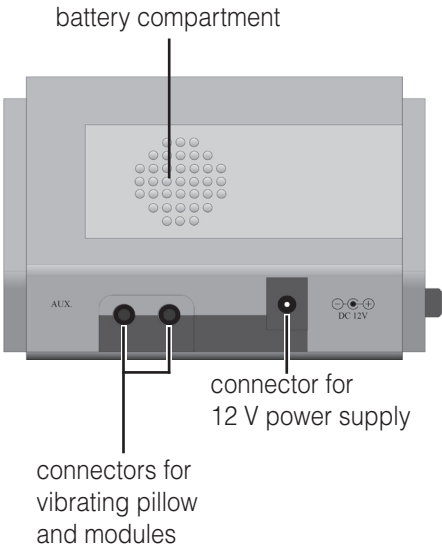
view of left side

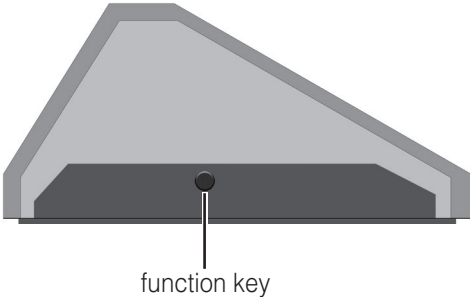Congratulations on the purchase of your »sign●lux« alarm monitor »alarmo«.

You have chosen a highly modern and reliable transmitter for the wireless »sign●lux« home signalling system.

Please read this manual thoroughly to be able to set up the unit correctly and to familiarise yourself with all of the features the system offers.

Package contents

Please check if all of the following components are included:

- »alarmo«
- rechargeable battery
- 12 V power supply unit
- manual
- warranty card

Should any parts be missing, please immediately contact your dealer or the manufacturer directly.

Function

»alarmo« has been conceived to enable the easy incorporation of already existing smoke detectors into the »sign●lux« signalling system, without need for modification or wiring.

»alarmo« scans environmental noises in proximity for sound patterns typical of smoke detectors and fire alarm systems of any type.

False alarms, possibly caused by alarm sounds not actually within in the room (e. g. passing emergency vehicles, TV and movie sound effects) are prevented by means of intelligent analysis and filtering of the registered sound.

If an alarm sound is identified as such, »alarmo« transmits a radio signal to all accordingly registered »sign●lux« receivers within range, which will indicate the incident with light, vibration or sound - depending on the type of receiver.

Signals and controls

Battery status indicator

A red illuminated battery status indicator means that the inserted rechargeable battery is almost depleted. The battery needs to be recharged.

A permanent green light indicates a charged battery and an operational »alarmo«.

Transmission status / momentary indicator

The transmission status indicator lights up when a radio signal is sent to the »sign●lux«-receivers.

Function key

The operational state of »alarmo« can be tested by pressing the function key.

After registering »alarmo« with at least one »sign●lux«-receiver (see page 10), the function can be tested by pressing this key and checking the alert given by the receivers.

Note (special accessories)

»alarmo« is able to trigger various indicator modules, e. g. you can have a light flash directly released by the »alarmo«.

The following modules (not included) can be plugged into to the AUX output connectors:

- Vibrating pillow (A-3303-0)
- Flash module MF-1 (A-3320-0)
- Switch module MS-1* (A-3330-0)
- Acoustic module MA-1 (A-3340-0)

*Only for use with Type C (CEE 7/16 Europlug) and Type F (CEE 7/4 „Schuko“) outlet connectors

Setting up the system

1. Placement

»alarmo« scans and analyses the noises within its environment for sounds typical of smoke detectors and similar alarm systems. In order to assist this function, »alarmo« should be placed as close as possible to the sound source of the alarm, which is to be monitored (e. g. on a shelf underneath a smoke detector).

In order to gain the best possible transmission range, please note the following recommendations:

- The transmitter should not be placed close to devices emitting intense interference, like TV screens, computer monitors or mobile phones.
- The transmitter should not be placed on objects made of metal, as the transmission range is diminished by the metallic shielding effect.

Note:

The device is not able of identifying, recording or even transmitting voices and spoken conversation.

Therefore there is no possible invasion or endangerment of the user's privacy.

2. Mains operation

Insert the plug at the end of the 12 V power supply into the corresponding connector on the rear side of the »alarmo«.

Then plug the power supply unit into a mains power socket in your room.

Power failure backup

Please note that the device is designed for mains operation.

The rechargeable battery is used as backup only in the event of a power failure.

The stand-by time of the battery is approx. 2 hours.

3. Inserting the rechargeable battery (included)

Open the battery compartment by sliding the cover to the right of the device using your thumb.

To connect the rechargeable battery, insert the battery plug into the outlet in the battery compartment.

The battery is charged when »alarmo« is connected with the mains power supply.

The battery has a charging time of 12 hours. Use and charge only the special battery A-2995-0.

4. Registering »alarmo« with the receivers

»alarmo« has to be registered with every »sign●lux«-receiver, which is later supposed to react to it. By doing so, the type of signalization is also defined.

To register the transmitter, press and hold the according receiver's function key.

The receiver's event symbols light up in sequence.

Release the function key as soon as the desired symbol is lit up. Now press the function key of the »alarmo«.

Wait for approx. 20 seconds and test the setup by pressing the »alarmo«'s function key. The receiver should now react, showing the according symbol.

Up to 6 different transmitters can be assigned to a single symbol.

Important notice:

»alarmo«'s radio alarm signal is triggered by the computer-assisted positive analysis of environment sounds.

Despite of digital filtering and programmed threshold values, the possibility of a false alarm cannot be fully excluded, due to the almost limitless variety of possible sounds (e. g. siren-like custom ringtones on mobile phones).

Due to these facts, the Humantechnik GmbH will take no responsibility for false alarms and their effects.

For this reason we advise to test the function of the »alarmo« by releasing a test alert on the alarm device (smoke detector, alarm system) which is to be monitored.

Sensitivity

»alarmo« has been developed and tested by considering and utilising many different types of alarm tones.

It is possible to alter the thresholds of your »alarmo«'s acoustic analysis, should it fail to respond to the monitored alarm device, even after repeated tests and alternating positions.

This involves opening the chassis and a certain process.

In this case, please contact the Humantechnik GmbH for advice.

Radio range

Note the following points when using the »alarmo«:

The radio signals are designed to have a range of 200 meters under perfect conditions. The following conditions may result in a reduced range:

- Signal transmission through building structures or vegetation.
- Interference from TVs, computer monitors and mobile phones. Under certain circumstances, this may prevent the receiver from functioning at all.
- Positioning of the transmitter at an unfavourable distance from surfaces which reflect radio waves, for example floors and walls. This may lead to the signal wave being weakened or even cancelled out by the reflected wave.
- Metallic objects reduce the range due to the shielding effect of the metal.
- Particularly in towns and cities, there may be many other sources of radio waves which can disrupt the original signal.
- Having units operating at similar frequencies located close to one another may also result in mutual interference.

Maintenance and care

»alarmo« does not require any maintenance. If the unit does become dirty, simply wipe it clean with a soft, damp cloth.

Note: Disconnect »alarmo« from the power supply!

Never use spirits, thinners or other organic solvents. Do not set »alarmo« up where it will be exposed to full sunlight for long periods. In addition, it must be protected against excessive heat, moisture and severe mechanical shocks.

Note: This product is **not** protected against splash water. Do not place any containers filled with water, such as flower vases, or anything with an open flame, such as a lit candle, above or near the product.

Please make sure that the batteries are not exposed to excessive heat, such as sunlight, fire or anything similar.

Warranty

The »alarmo« is a very reliable product. Should a malfunction occur despite the unit having been set up and operated correctly, please contact your dealer or the manufacturer directly.

This warranty covers the repair of the product and returning it to you free of charge. It is essential that you send in the product in its original packaging, so do not throw the packaging away. The warranty does not apply to damage caused by incorrect handling or attempts to repair the unit by unauthorised persons (destruction of the seal on the unit).

Repairs will only be carried out under warranty if the completed warranty card is returned accompanied by a copy of the dealer's invoice/till receipt.

Always specify the product number in any event.

Disposal of used electric and electronic units (applicable in the countries of the European Union and other European countries with a separate collection system).

The symbol on the product or the packaging indicates that this product is not to be handled as ordinary household waste but has to be returned to a collecting point for the recycling of electric and electronic units. You protect the environment and health of your fellow men by the correct disposal of this products. Environment and health are endangered by a faulty disposal. Material recycling helps to reduce the consumption of raw material. You will receive further information on the recycling of this product from your local community, your communal disposal company or your local dealer.

Technical data

Power supply:	12 V-power supply unit at 110 - 230 V, 50-60 Hz or via rechargeable battery A-2995-0
Rechargeable battery standby time:	Approx. 2 hours
Frequency:	868.35 MHz (915 MHz in according export countries)
Operating range:	0° C to 40° C
Transmission indicator:	1 red LED
Battery level indicator:	1 red/green LED
Height:	65 mm
Width:	100 mm
Depth:	135 mm
Weight:	220 g

This device satisfies the following EU directives:


- 2011/65/EG RoHS directive
- 2002/96/EG WEEE directive
- 1999/5/EG R&TTE directive


Compliance with the directives listed above is confirmed by the CE seal on the device. CE compliance declarations are available on the Internet at www.humanteknik.com.

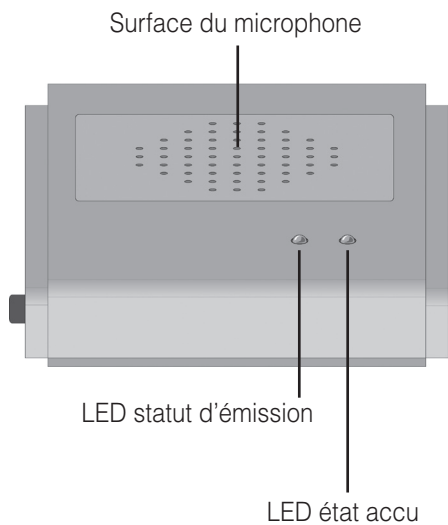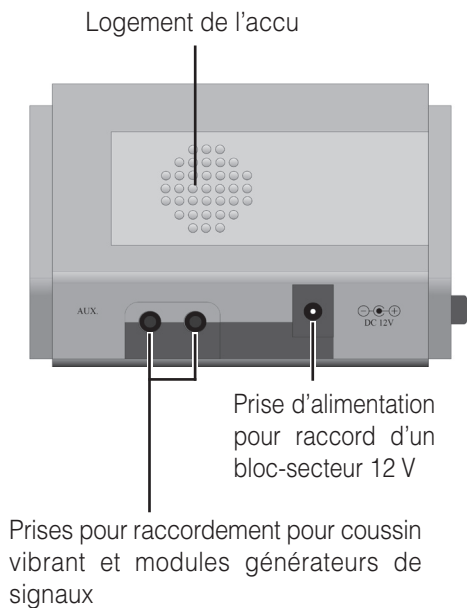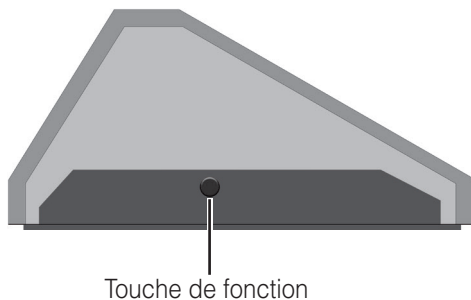
Technical specifications subject to change without prior notice.

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Sarabec Ltd